

The LiveLife Mobile Personal Alarm.



Perfect for people at risk of falls

Texts up to 6 emergency contacts with your GPS location and then calls them consecutively. Can call 000.

Speak hands free through your pendant.

After activation the pendant makes audio announcements to guide the wearer through their emergency.



Finally an emergency alarm that works anywhere.



Bluetooth 5

How it works.

Revolutionary technology



The **LiveLife 4Gx Mobile Alarm** is an SOS alarm that works wherever there is 4Gx, 4G or 3G coverage through the Telstra mobile network which covers 99.5% of Australia's population.



When the help button on the pendant is pressed it uses Bluetooth and GPS to send out up to 6 help text messages with the wearer's location on **Google Maps**. This is accurate to 2 metres. When family or friends receive the call from the Mobile Alarm they will already know exactly where the wearer is located. After activation the pendant makes audio announcements to guide the wearer and keep them informed.



The Mobile Alarm then begins the call sequence. When the first contact answers the pendant allows the wearer to speak and listen hands free. It's just like a normal phone conversation. The pendant has loud and clear audio. When the wearer is finished talking they can end the call very easily.

Other contacts who have received the help text message can then try calling the pendant if they like and speak to the wearer.

At any time the contacts can find out the location of the pendant by sending a simple text to it. It will reply with its exact location on Google Maps.



www.livelifearms.com.au



Whether you're at home or outdoors and need help, simply press the SOS button.



It then sends 6 help texts with your exact location via GPS and then calls those people consecutively.



When one of your emergency contacts answers you can talk to them through your pendant.

LiveLife Mobile Alarm


Safety and independence for active seniors.

With the **LiveLife** 4Gx Mobile Alarm you can go anywhere knowing help is just a press of a button away. It is perfect for people who want to feel safe everywhere, not just at home.

It is a stylish pendant that works on the Telstra Mobile network but without the complicated features of mobile phones. The wearer doesn't need a mobile phone, the NBN or the internet.

When activated, the pendant sends up to 6 help messages via text using the built-in Swiss precision GPS with a link to **Google Maps** showing the wearer's location to within 2 metres.

It then starts calling up to 6 emergency contacts one after the other using the Telstra Mobile network. You can even include 000 as a contact!

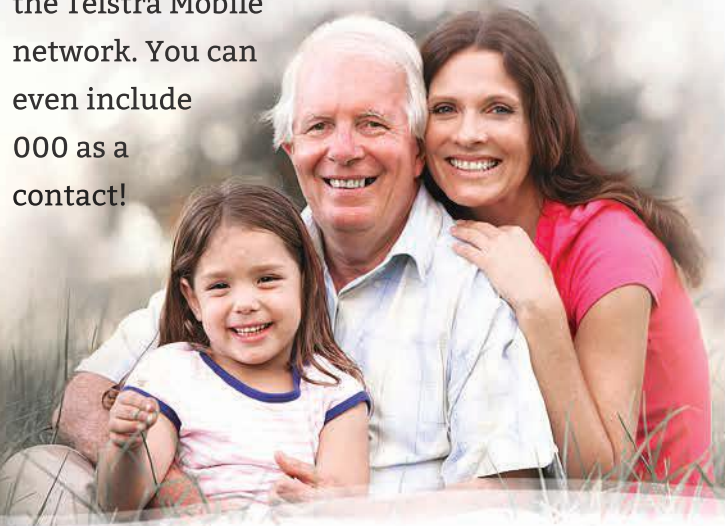

When one of the emergency contacts answers, the wearer is able to talk hands free through the showerproof pendant. Anyone who knows the mobile number of the pendant can call it and it answers in speaker-phone mode automatically.

You can easily locate the pendant via text and it will send you the location on Google Maps.

 **Perfect for people** who are at risk of falls.

Includes a **fall detection** capability designed to activate for more serious falls likely to result in unconsciousness. It alerts all contacts.

No monthly fees ever.
One off cost only.

In an emergency it contacts family and friends, not strangers in a call centre who will just call your family and friends anyway.

Ready to use right out of the box. No setup needed.

It has its own SIM card installed and mobile number. It comes all setup, ready to wear with 12 months warranty and a 14 day moneyback guarantee.

Small, lightweight and comfortable to wear.

The LiveLife 4Gx Mobile Alarm pendant weighs the same as 3 twenty cent pieces. It has a soft, non-slip surface and comes with a 'safety' lanyard, belt clip and wristband.

Being showerproof you can wear it in the shower or even a tropical downpour.

The pendant fully recharges in the charging station in 25 minutes.

The wearer can also call one of the emergency contacts at any time if they wish.



000
CAN CALL
EMERGENCY
SERVICES

GERMAN
DESIGN &
ENGINEERING

SWISS
PRECISION
GPS

reddot award
product design
winner



Intelligent. So simple to use.

Smart charging station

The **LiveLife Mobile Alarm** comes with a new Smart Charging Station. It's Bluetooth automatically stays connected to your pendant when you are at home so you can be sure the alarm knows where the wearer is at all times.

The Smart Charging Station has it's own SOS button and loudspeaker. The Station talks to you just like the pendant does to tell you when the device is charging and when your pendant is activated it guides you through your emergency.

With the new Bluetooth 5 technology the battery life of the pendant is extended to four to five days of use. It only takes 45 minutes to charge in the charging station.

The Smart Charging Station has it's own battery inside. If there is a power blackout and you alarm's battery is low you can still charge it.

Customised set up

For all our customers we custom program the alarm which includes recording the wearer's residential address into the Smart Charging Station for location accuracy in an emergency.

 **Bluetooth 5**



Smart Talk Feature



How to order

Order online.



You can order securely online at our website:

www.livelifearms.com.au

Order by phone.



Freecall: 1800 936 774

More information

Visit our website.



www.livelifearms.com.au

info@livelifearms.com.au

Cost

Ready to go.



Just one affordable payment to own outright* No monthly monitoring centre fees. Programmed and **ready to use straight out of the box.**

Device uses a SIM card on the Telstra network. First year of unlimited emergency use is included in purchase price.

*45 per year in 2nd year and each year following. No contracts.

See website for full details.



Freecall: 1800 936 774

Web: www.livelifearms.com.au

Email: info@livelifearms.com.au

Address: 4/7 Pemell Street (Above Woolies), Toronto NSW Australia 2283